



# MMT Bladder Scanner




Real-time, save, easy  
and non-invasive bladder  
scanning.

- One button controls the scanner start and stop.
- Manual bladder correction mode.
- Physiological anomalies that may affect volume calculation are visualized.

MMT Bladder Scanner software is easy  
to operate

- 3D volume calculations – PVR saved into Uromic application automatically.
- PVR is available in seconds.

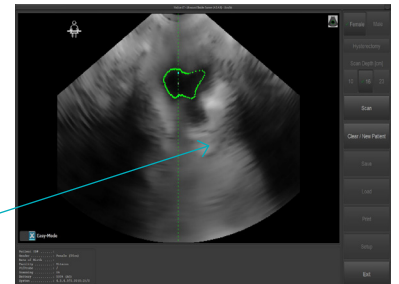

## Special holders for Uromic systems



Uromic QuickStep

Uromic Samba

Uromic Jive

Uromic Jive



# MMT Bladder Scanner Stand alone system

MMT Bladder Scanner measures ultrasonic reflections within a patient's body and differentiates the urinary bladder from the surrounding tissues.

- No sonographer skills needed.
- Fully automatic multi-plane scans from one position on the body.
- Auto location and auto tracking of the bladder in real-time.
- Advanced signal and image processing for sharp and clear images.
- Full-page reports in JPEG and Dicom.
- Post-test quality indicators help to train users for optimum results.

## THREE SCAN DEPTHS

- |         |                       |
|---------|-----------------------|
| • 10 cm | designed for Children |
| • 16 cm | designed for Adults   |
| • 23 cm | designed for Obese    |

## TECHNICAL SPECIFICATION

- |                   |                               |
|-------------------|-------------------------------|
| • Volume range    | 0 – 1000 ml                   |
| • Accuracy        | +/- 10% of reading, +/- 20 ml |
| • Scanning method | Sector, 180 degrees           |

Through USB connector it's possible to connect the MMT Bladder Scanner to any tablet, laptop or PC.